

# OrthoLavage<sup>®</sup>

## PULSE LAVAGE SYSTEM

Removable splash shield.

Spray stem and suction tube can be removed independently for greater versatility.

Fracture-resistant spray stem and suction tip.

High pressure lavage system.

The OrthoLavage<sup>®</sup> system provides:

- Superior cleaning power
- High pulse frequency
- Increased pulse strength
- Optimal saline flow rate

Fully variable speed trigger for optimal control of flow and pressure.

Ergonomic hand grip.

No batteries, gives more power and less waste


Integrated irrigation and suction tubing.  
Non-clogging suction design.

Single use Lavage System



# OrthoLavage®

PULSE LAVAGE SYSTEM

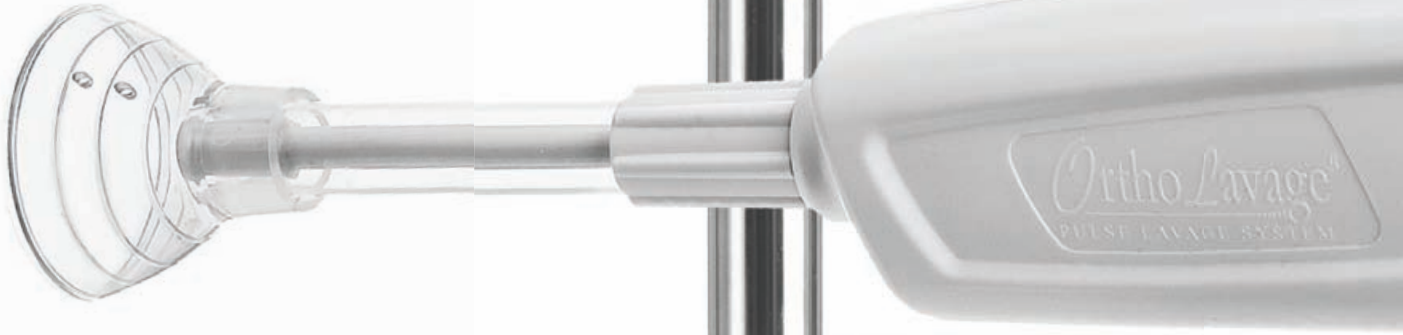

## *Superior Suction Design*

Excellent, non-clogging suction design ensures smooth passage of debris.


## *No Batteries - More Power*

OrthoLavage® Power Supply with three mounting options;

- Table Mount
- I.V. Pole back mount
- I.V. Pole hang mount

LED indicators for function monitoring.



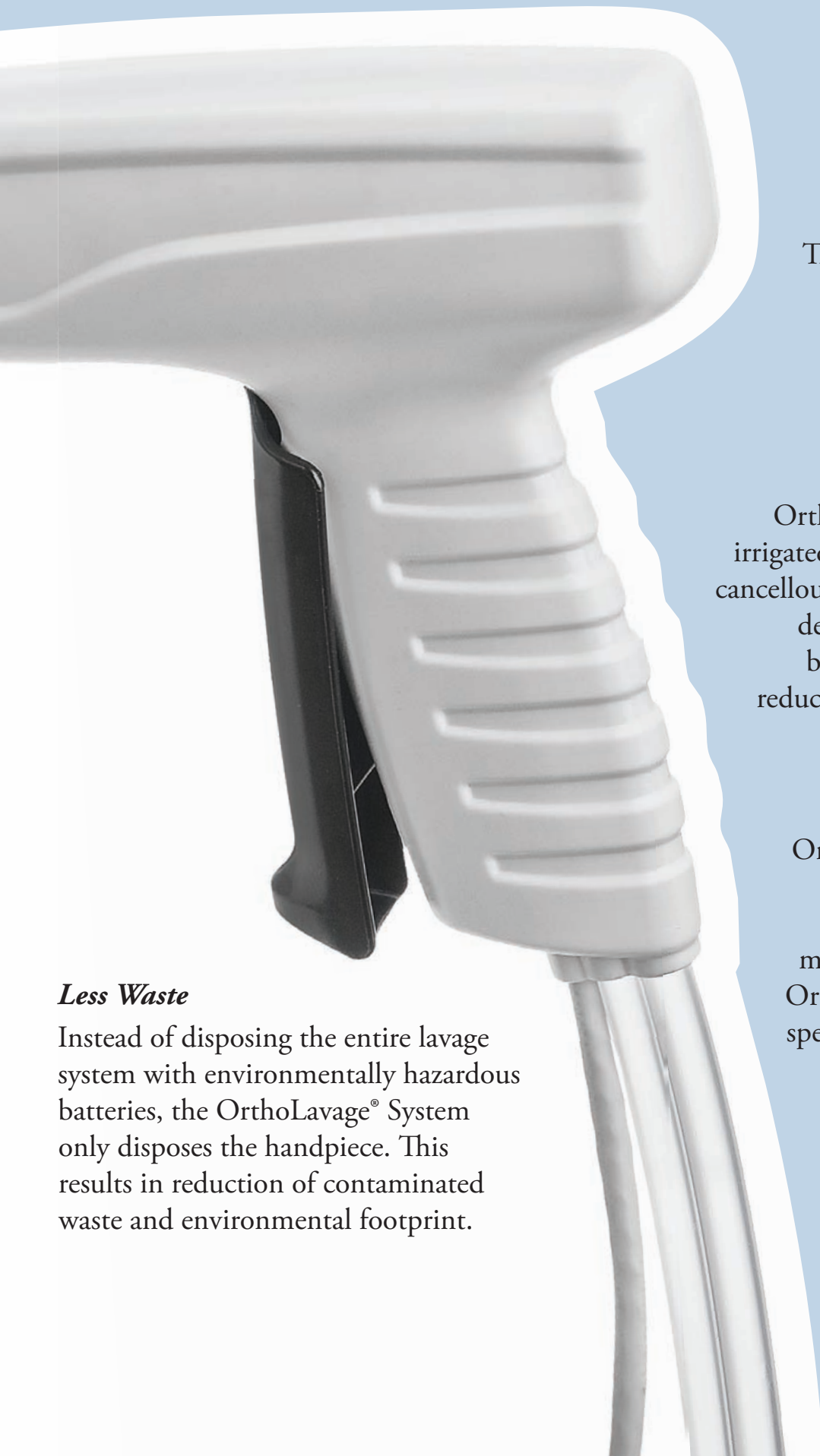

## *Technology*

The OrthoLavage® equipment is used both in Orthopedic surgeries and in Wound Management.

### *Orthopedic Surgery*

In Orthopedic surgery the OrthoLavage® uses pulsed jets of irrigated saline to penetrate into the cancellous bone, removing blood and debris, supplying an improved bone-to-cement interface and reducing the risk of fat embolism.

### *Wound Management*

In Wound Management the OrthoLavage® facilitates quick, efficient removal of necrotic tissue, bacteria and foreign material from the wound. The OrthoLavage® Soft Tissue Set is specifically designed for wound management providing a gentle pressure, optimal for cleaning the wound with minimal harm to healthy tissue.

### *Less Waste*

Instead of disposing the entire lavage system with environmentally hazardous batteries, the OrthoLavage® System only disposes the handpiece. This results in reduction of contaminated waste and environmental footprint.

# Single use Lavage System

## Disposables

---

OL1100

OrthoLavage® Hip/Knee Set

OL1200

OrthoLavage® Soft Tissue Set

## Power Source

---

OLAC2230

OrthoLavage® AC Power Supply 230 V

OLAC2115

OrthoLavage® AC Power Supply 115 V



Manufactured for:  
Pulse Lavage AB  
Stabby Allé 2  
SE-75229 Uppsala  
SWEDEN  
E-mail: [pulselavage@pulselavage.se](mailto:pulselavage@pulselavage.se)  
Phone: +46-18-55 55 05  
Fax: +46-18-51 50 50  
Web: [www.pulselavage.se](http://www.pulselavage.se)

Manufactured by:  
Oy Fluorplast AB  
Västervägen 8  
66240 Petalax  
FINLAND

  
0537

RCD No  
001934589